VOCO

# VOCO Ionofil® Molar

GB

CE 0482

In accordance with DIN EN ISO 9917-1

## Instructions for use

### Product description:

**VOCO Ionofil Molar** is a radiopaque glass ionomer cement for fillings with chemical adhesion to enamel and dentine and continuously high release of fluoride ions.

**VOCO Ionofil Molar** is available in the shades A1, A2 and A3.

### Indications:

- Restoration of V-shaped defects and cervical enamel or root erosions including class V cavities
- Restoration of deciduous teeth (especially class I)
- Fillings of class II cavities (not occlusion-bearing) by means of tunnel preparation
- Build-up fillings and linings, core build-up
- Temporary long-term treatment of class I and II cavities

### Method of application:

#### Preparation:

The adhesion of **VOCO Ionofil Molar** to enamel and dentine prevents the formation of marginal gaps and creates retentions. Undercuts give additional hold to the filling. The layer thickness of fillings should not be less than 1 mm. Avoid feather-shaped margins. In case of doubt choose conventional preparation. Clean the cavity (e. g. 3 % hydrogen peroxide), thoroughly rinse with water and dry (warm air). Do not overdry dentine. If the remaining dentine is thicker than 1.5 mm additional lining is not necessary. Cover areas in close proximity to the pulp with calcium hydroxide to avoid pulp irritations which may occur in cases of insufficient dentine thickness.

Apply **VOCO Ionofil Molar** at room temperature (19 °C - 23 °C); in case of high temperatures the use of a cooled glass plate is advantageous. Low temperatures prolong the working time.

Furthermore, **VOCO Ionofil Molar** is suitable for use for the purposes of ART-technique (Atraumatic restorative treatment).

#### Mixing powder/liquid:


Shake powder well before use!

Mix one level measuring scoop of powder with one drop of liquid. This corresponds to a mixing ratio of 4.0 parts powder to 1 part liquid by weight. Hold dropping bottle in vertical position and place 1 drop of liquid onto a dry glass or plastic plate. Hold bottle high enough so that a single drop can fall down freely. Only use drops without enclosed air bubbles!

Mix the powder into the liquid in two equal portions using a clean metal spatula, or preferably with a clean **plastic spatula**. Only mix in the second portion when the first portion has been mixed in homogeneously. The mixed cement has a firm, packable consistency.

Avoid any contact with humidity!

Never add new liquid during mixing!

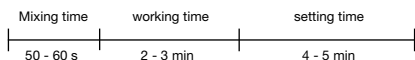The mixing time takes approx. 50 - 60 s and should not exceed 1 min.

#### Application:

After application into the cavity **VOCO Ionofil Molar** will set within approx. 2 - 3 min. During the placing of a filling the working area must be kept dry. The use of a rubber dam is recommended. Linings or build-up fillings can be carefully finished after setting (approx. 4 - 6 min after application).

Glass ionomer cements are susceptible to humidity (breath, saliva) or, respectively, dehydration during the setting time, especially directly after application. If possible, cover immediately with a strip or matrix or apply a protective varnish.

Finishing and polishing (flexible polishing discs) of the fillings should follow 7 min after application at the earliest. Afterwards, cover again with varnish. The best finishing results will be achieved approx. 24 h after the filling has been inserted.

#### Information, precautionary measures:

Close powder tightly after use in order to keep it dry and to maintain its high physical strength. Remove any remaining liquid from the outside of the liquid dispenser (by tissue) to keep opening clean for a correct mixing ratio. **VOCO Ionofil Molar** should only be mixed immediately before use to ensure that the entire working time is available for placing and modelling. The development of fissures due to delayed use of the mixed cement can thus be prevented.

**VOCO Ionofil Molar** contains polyacrylic acid, fluorosilicate and parabenes. **VOCO Ionofil Molar** should not be used in case of known hypersensitivities (allergies) to any of these ingredients.

#### Storage instructions:

Store **VOCO Ionofil Molar** powder and liquid containers tightly closed at room temperature (15 °C - 23 °C) and protected from light. Do not use after expiry date.

Our preparations have been developed for use in dentistry. As far as the application of the products delivered by us is concerned, our verbal and/or written information has been given to the best of our knowledge and without obligation. Our information and/or advice do not relieve you from examining the materials delivered by us as to their suitability for the intended purposes of application. As the application of our preparations is beyond our control, the user is fully responsible for the application. Of course, we guarantee the quality of our preparations in accordance with the existing standards and corresponding to the conditions as stipulated in our general terms of sale and delivery.

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